

NEPTUNE

MULTIFUNCTIONAL POOL LIFT



nibotechnics

The Neptune pool lift gives your patient a safe, comfortable and dignified way of transferring from the poolside into the water. It is also easy and effortless for you, the operator, to use.

The system includes chairs and comfortable stretchers which attaches easily to the lift. The Neptune pool lift suits above ground pools as well as deck level pools.


TRANSFER AND LIFTING SYSTEM FOR SITTING PATIENTS



The Neptune chair transports the patient and attaches to the lift. The fold-back armrests give him additional security. After the chair is attached to the lift, remove the sub-chassis. The swivelling castors will help you to do this smoothly.

ALTERNATIVE CHAIRS

In addition to the bucket seat, a seat with a commode aperture is available. The commode facility also promotes easy washing in the shower when the chair is provided with the additional sub-chassis. In addition, there is also a swimchair, suitable to go with the fixed sub-chassis in the water.

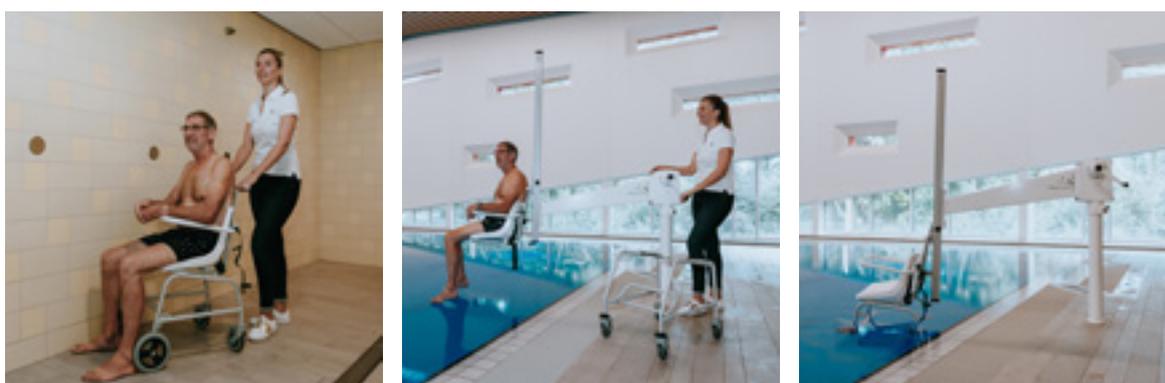

All chairs may be used in the changing room, in the shower and the commode chair for toileting too.

Using more than one chair makes the transportation of several patients to the pool quick and efficient.


The swimchair has a permanently attached sub-chassis which makes it possible for your patient to remain seated on the chair for exercises.

TRANSFER AND LIFTING SYSTEM FOR RECUMBENT PATIENTS


There are two stretchers available for recumbent patients to bathe; one with a comfortable stretcher platform and quilt and one with a water-permeable fabric. Both stretchers are available with a fixed or a height adjustable trolley.



The stretcher is wide and comfortable, and the trolley has swivelling castors which make it easy to operate.

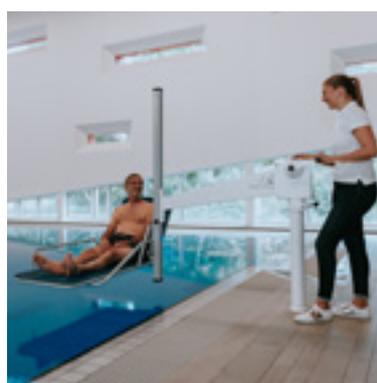


Attach the stretcher to the lift. Detach the trolley. Lower your patient into the water. This is a smooth, safe operation.

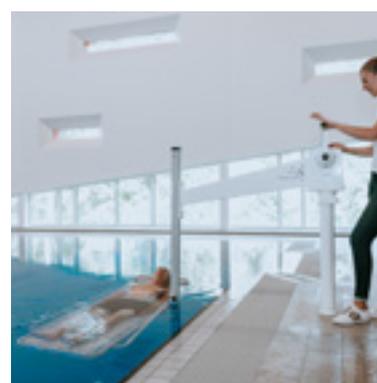


Your patient can stay on the stretcher for his exercises, or, with some assistance, ease himself into the water.


ACCESSORIES

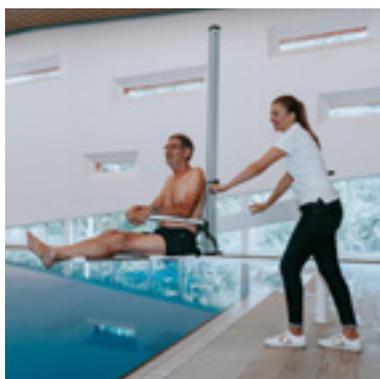

A legrest for the bucket and commode chair makes the transfers into above ground pools smoother.

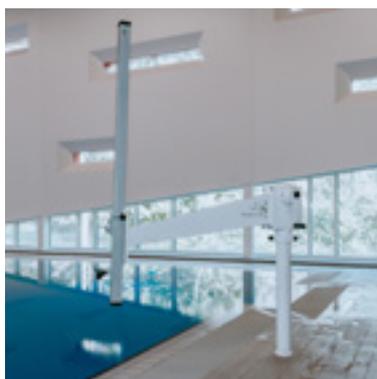

The Neptune pool lift can be secured to the floor either by means of a socket or by a base plate.

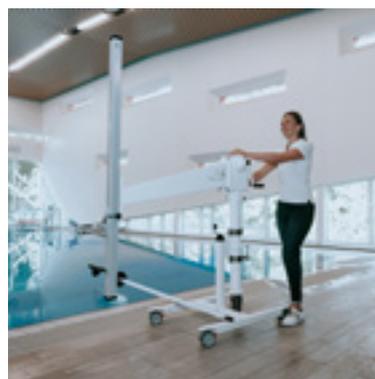

The transporter makes lifting and moving the Neptune to a store, or from one installation point to another in the pool area, virtually effortless.



Chairs are available with armrests which fold back. Safe at the transfer and comfortable when getting on and off the chair.



The removable chassis is equipped with brakes on all four castors to make the transfer to and from it safe and simple. It is available with a foot rest or a foot bar.



The fixed chassis of the swimchair is made of high quality materials for use in water. The large front wheels make it easy to move in the pool.

PRODUCT SPECIFICATIONS

safe working load	160 kg (352 lbs)
lifting stroke	1270 mm (50")
lifting stroke above water surface, min.	680 mm (26 3/4")
lifting stroke above water surface, max.	590 mm (23 1/4")
total product height	2340 mm (92")

Nibotechnics nv

Heikant 5 | 3930 Hamont-Achel | Belgium

T +32 (11) 80 20 40 | info@nibotechnics.com | www.nibotechnics.com