

EMERGENCY CARTS

Rapid and organized response to a cardiac emergency is essential. The Waterloo aluminum cart is about 30% lighter than a standard steel cart, offering faster response times. The included steering caster allows the cart to take tight corners further improving maneuverability.

Work surface and drawer organization are vital during a Code. Our dual slide out shelves virtually double your cart top work surface, allowing you to move equipment/supplies to the cart top and within easy reach.

Inside our drawers you'll find durable metal drawer dividers that will keep supplies organized and orderly. No more fumbling around for items in the heat of an emergency. Our metal dividers will stand up to years of use and can be sub divided with optional plastic sub dividers.

Dual slide out work surfaces operate independently and virtually double your work surface!

Included plastic top with cutaway front is designed to be "Wrist friendly" as you work on the cart top.

Large stylish/contemporary drawer handles offer a more updated look and enhanced functionality

Adjustable metal dividers (One long and One short) included in all aluminum cart drawers.

Advanced nylon reinforced polycarbonate bumper provides a wider wheel base and greater cart stability and protection. Single piece construction resists scuffs and scratches.


5" Premium "Easy roll" sealed bearing casters with precision bearing swivels. Thread guards minimize sutures and other materials from collecting around axles and keep carts rolling smoothly.



UTRLA-333369-RED

RECOMMENDED

IV POLE
Standard style IV Pole.
Part# IV-2
SEE PAGE 30.

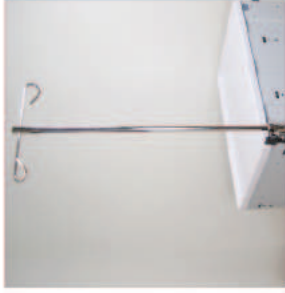

OXYGEN TANK BRACKET
Powder coated oxygen tank bracket.
Part# OTB-5
SEE PAGE 30.



CARDIAC BOARD
Slim cardiac board and mounting brackets.
Part# CBB-9
SEE PAGE 30.



FEATURE

ELITE ALUMINUM

- Unicart Construction
- Ball Bearing Slides
- Easy-roll Casters
- Includes Drawer Divider Set
- Universal Accessory Mounting System
- Dual Pullout Work Shelves
- Lightweight Construction

WRAP AROUND RAILING

Cart top railing keeps items on cart top during transport.

Part# R-4 (aluminum carts)

R-3 (Steel carts)

R-JR (Med-Jr carts)

SEE PAGE 32.



SHARPS CONTAINER BRACKET

Accommodates standard 5QT sharps container. Locking version also available.

Part# SB-2

SEE PAGE 30.



SWIVEL MONITOR SHELF

Monitor/Defibrillator shelf swivels 360 degrees and locks into position to accommodate defibrillator access from any side of the cart.

Part#

ER-SHLEA (Aluminum Carts)

ER-SHLEF (Steel carts)

ER-SHLE2000 (ER2000 carts)

SEE PAGE 32.



COMPARISON

TRADITIONAL STEEL

- Unicart Construction
- Ball Bearing Slides
- Easy-roll Casters
- Includes Drawer Divider Set
- Universal Accessory Mounting System
- Dual Pullout Work Shelves
- Lightweight Construction

EMERGENCY CARTS

Economical steel crash carts offer a robust work platform that can be configured with accessories to address your code cart needs.

Collapsible side shelves are available for added work surface and a variety of drawer trays and dividers are offered to further organize your cart.

Shares many of the same great features as our aluminum cart but at a more economical initial price.

If you are considering a steel cart and are looking to organize it with plastic trays/dividers as well as a side shelf or two, you may wish to consider our aluminum cart as the included dividers and side shelves can offset the added cost of upgrading to aluminum.



Smooth, easy-to-clean and removable work surface.

Pre-drilled holes to mount accessories.

Advanced nylon reinforced polycarbonate bumper offers a wider wheel base for more stability while offering superior protection. Single piece construction that resists scuffs and scratches.

5" Premium "Easy roll" sealed bearing casters with precision bearing swivels. Thread guards minimize sutures and other materials from collecting around axles and keep carts rolling smoothly.

UTRLU-433369-RED

CRASH CARTS


Crash carts come in all shapes and sizes, and at Waterloo we acknowledge that fact by manufacturing a wide variety of cart sizes designed to fit your specific needs. Simply choose your desired drawer combination of 3", 6", 9" or 12" drawers and we can help you pick the correct cart height to accommodate your configuration. From there, you can choose between an all red or all blue cart depending on your facility's designated crash cart color.

If you'd like to color code drawers to make them easy to identify quickly (example: White airway drawer), we are happy to accommodate that request at no extra charge! Just tell us what drawer front and what color and we can

personalize your carts for you. Another advantage of working with a true manufacturer of medical carts!

We can also outfit your crash carts with individual breakaway gatelocks over various drawers. These can not only save re-stocking time, but add extra security to drawers that need it (See page 6 for lock options).

Waterloo aluminum crash carts available in Extra-Tall, Tall, Medium and Short cart sizes (See page 5 for sizing chart).

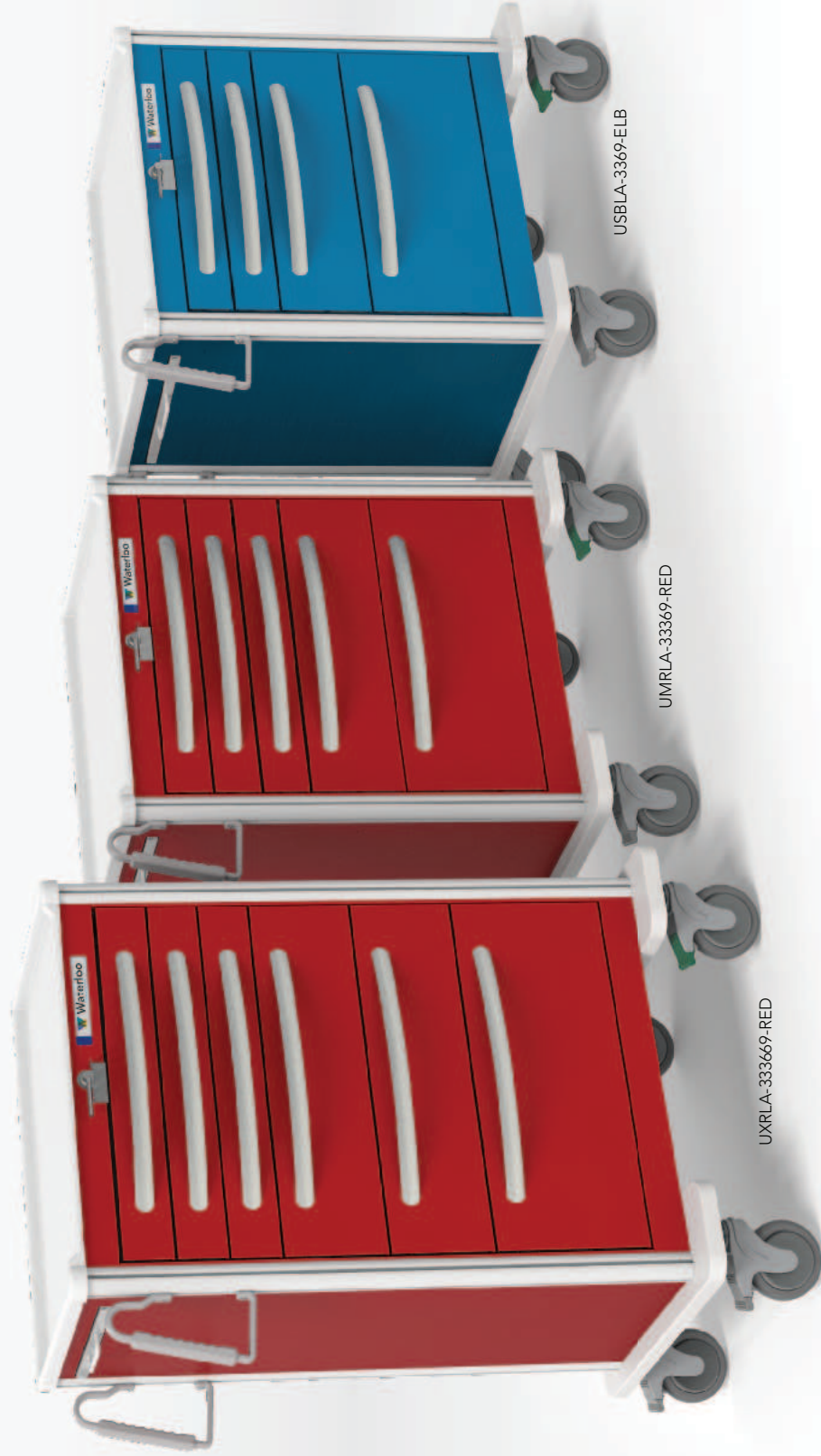

CRASH CARTS

Steel carts come in many of the same cart heights as our versatile aluminum carts. Just like our aluminum carts, choose your desired drawer combination of 3", 6", 9" or 12" drawers which will help us determine the correct cart height for your configuration. Then, choose between an all red or all blue cart depending on your facility's designated crash cart color.

Should you like to color-code any drawers, simply let us know which drawer and what color (from one of our 16 available colors) and we are happy to accommodate your request. We never charge for this option!

Individual gatelocks are available on steel carts and can be installed at our factory for you.

Waterloo steel crash carts available in Extra-Tall, Tall, Medium and Short cart sizes (See page 5 for sizing chart).

The Waterloo ER-2000 takes an innovative approach to crash carts and is designed to offer full access to all drawers in a code with the push of a button. Call for more details!



One press of the centrally located latch opens the dual front doors exposing all drawers for rapid access.

Drawer trays can be removed to be positioned where needed during a code.