

Functions and features	
White balance	Balances the color temperature in your image.
HbE Enhancement	Highlights the mucosal capillary morphology
Image Enhancement	Structure/Edge/Color
IRIS	Average , Peak, Automatic
Storage	USB
Color adjustment	Adjustment of Y(Luminance), R(Red), or B(blue) component within 21 levels (from -10 to 10)
Freeze and playback	Freezes a real-time image and plays it back
Digital zoom	1/1.2/1.5/2x

HD-GT/CL								
	Model	Working length	Diameter of distal end	Diameter of insertion tube	Diameter of instrument channel	Depth	View of field	Range of bending
Video Gastrosopes	HD-GT1	1050mm	9.6mm	9.6mm	2.8mm	2-100mm	140°	U210°D90°L100°R100°
	HD-GT1P	1050mm	9.0mm	9.0mm	2.8mm	2-100mm	140°	U210°D90°L100°R100°
Video Colonscopes	HD-CL1I	1330mm	12.8mm	12.8mm	3.7mm	2-100mm	145°	U180°D180°L160°R160°
	HD-CL1L	1680mm	12.8mm	12.8mm	3.7mm	2-100mm	145°	U180°D180°L160°R160°

 **+49 (0) 8104-8908-0**

Headquarter Europe (Aohua affiliate company)
WISAP Medical Technology GmbH
Fichtenstrasse 27
85649 Brunnthal-Hofolding
Email: Info@aohua.eu

www.Aohua.eu

See the Manual for details of contraindications or precautions.

v1.0_2308



AC-1
Video Endoscopy System

Innovative compact system with above-standard configuration to explore a new horizon


One step ahead through innovation. The AC-1 is a compact system that combines several first-class technological advances. Experience the comfortable and price-efficient system and convince yourself of the high quality the AC-1 is offering.

Innovative solution



Compact design

In order to save the valuable but limited space on the cart for other peripherals, Aohua offers a compact design with a LCD touch panel on the front of the AC-1.

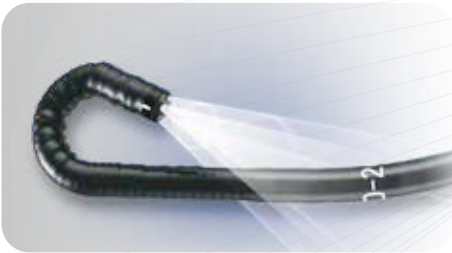

Distal end LED light source

The LED light source with a lifetime of 10,000 hours is integrated at the end of the distal tip, allowing to have a slimmer outer diameter. Even after long-term use, there is no need to compromise on brightness and consistency.



HD image quality and more

With the HD CMOS sensor, the AC-1 features true high-definition images for diagnostics. Combined with advanced digital processing technologies such as hemoglobin enhancement (HbE), digital chromoendoscopy.



Premium features available



Fully waterproof and one-touch connection

The waterproof one-touch connection brings convenience to your team and safety to your devices. Make use of the advantages the AC-1 is giving you.

